


Art.	Descrizione	Description
01CI4540	<p>PLASTICIMETRO DI PFEFFERKORN</p> <p>Strumento idoneo alla valutazione della lavorabilità di materiali ceramici. La misurazione avviene attraverso la caduta di un piatto calibrato, sul sottostante provino. Il provino viene formato con apposito attrezzo . Il plasticimetro ha due scale di lettura: una misura la deformazione in mm, ed una misura la deformazione del provino in accordo alle esperienze di Pfefferkorn. Fornito completo di attrezzo speciale per la formazione dei provini. Dimensioni di ingombro cm 26x38x54. Peso netto: kg 19</p> <p>Ricambi:</p>	<p>PFEFFERKORN PLASTICITY TESTER</p> <p>Instrument for evaluating the workability of ceramic materials. The measurement is based on a calibrated plate falling onto an underlying test body which is shaped by a special tool. The plasticity tester has two reading scales: one measures deformation in mm, while the other determines deformation of the test body according to Pfefferkorn tests. Supplied complete with special tool for shaping test pieces. Overall dimensions: 26x38x54 cm Net weight: 19 kg.</p> <p>Spare parts:</p>
01CI4540/A	<p>Attrezzo completo per la formazione del provino.</p>	<p>Complete set of equipment for shaping test pieces.</p>

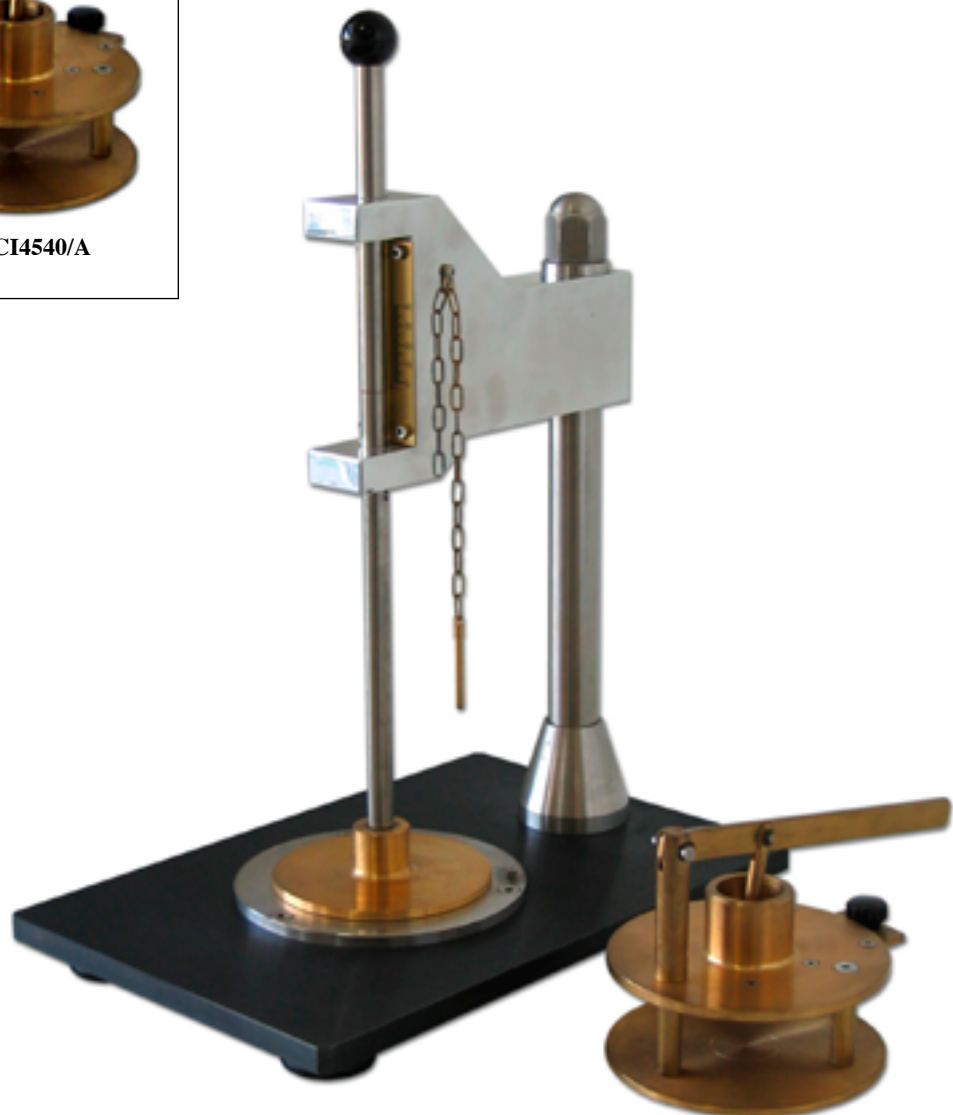

Art. 01CI4540